



PATIENT WARMING + GROUNDING

Warming Mattress + Return Electrode¹

Compare Technology to Megadyne™ Mega Soft™ Patient Return Electrode Pads²



More Info


U501

20 in x 32 in (50.8 cm x 82 cm)
2.0 lb (0.9 kg)



U502

20 in x 50 in (50.8 cm x 127 cm)
3.3 lb (1.5 kg)

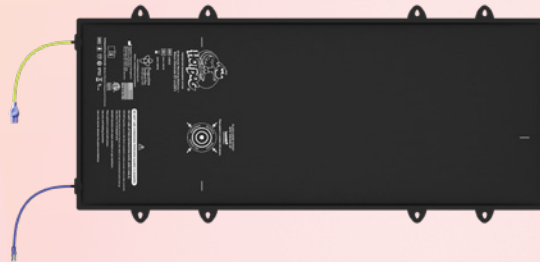

U530

20 in x 35 in (50.8 cm x 89.0 cm)
2.3 lb (1.0 kg)



Easy to Use

- Compatible with most electrosurgical generators³
- Use with one or two electrosurgical generators
- Lightweight for quick & easy set up and cleaning
- One product, two functions
- Requires no skin preparation, such as shaving
- Tattoos, piercings, and implants do not interfere with return electrode performance

Improved Safety

- Eliminates risk of burns from loose adhesive pads
- No adhesive - ideal for patients with delicate skin
- Capacitive-coupling return electrodes have been used for 20+ years in over 100 million procedures⁴

Patient Warming

- Air-free, safe, conductive heat
- Better temperature outcomes when combining a Mattress and Blanket

Cost Savings

- Reusable for up to 2 years
- Saves money and time compared to single-use pads or other reusable return electrodes
- Reduces SKUs and ordering hassles

Environmentally Friendly

- Reduces waste from single-use adhesive pads and single-use warming devices
- Lowers pollution from shipping and transport
- Meets green initiatives⁴

1. FDA cleared, 510(k) K230866, 510(k) K232627

2. Megadyne and Mega Soft are trademarks of J&J Ethicon

3. Adapter may be required

4. jnjmedtech.com/en-US/product/megadyne-mega-soft-patient-return-electrode-pad



Warming Mattress + Return Electrode Cables

A112

**Warming Mattress
Cable**


16ft (5 m)



A136

**Return Electrode
Cable**


13 ft (4 m)

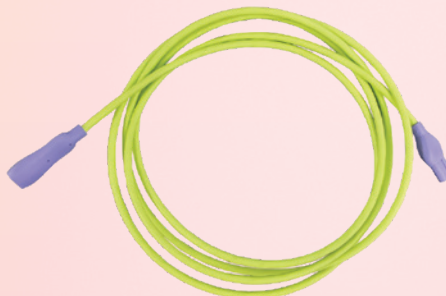

A137

**Return Electrode, Dual
Generator “Y” Cable**

102 in (2.6 m) (Each Y Length)
Total Length: 165 in (4.2 m)

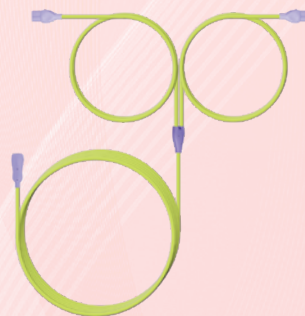

A138

**REM Enabled Return
Electrode Adapter**

16.5 in (42.0 cm)



A139

**Return Electrode
Adapter Cable**

16.5 in (42.0 cm)



**Augustine
Surgical Inc.**
Making Surgery Safer

www.hotdogwarming.com

©2026 Augustine Surgical, Inc. All Rights Reserved.

Patents: hotdogwarming.com/patents

Rx Only

M137EN Rev C 7/26

6581 City West Parkway
Eden Prairie, MN 55344
952-465-3500 or 866-484-3505
email: cs@aug surg.com