



WaffleGrip® Trendelenburg Positioning System + BackSaver™ SlideSheet



WaffleGrip
Warm & Secure




Learn More

Technology Integration


That Optimizes Surgical Efficiency and Safety

- Safe Trendelenburg Securement
- Efficient Patient Warming
- Easy Patient Positioning
- Short Turnover Times

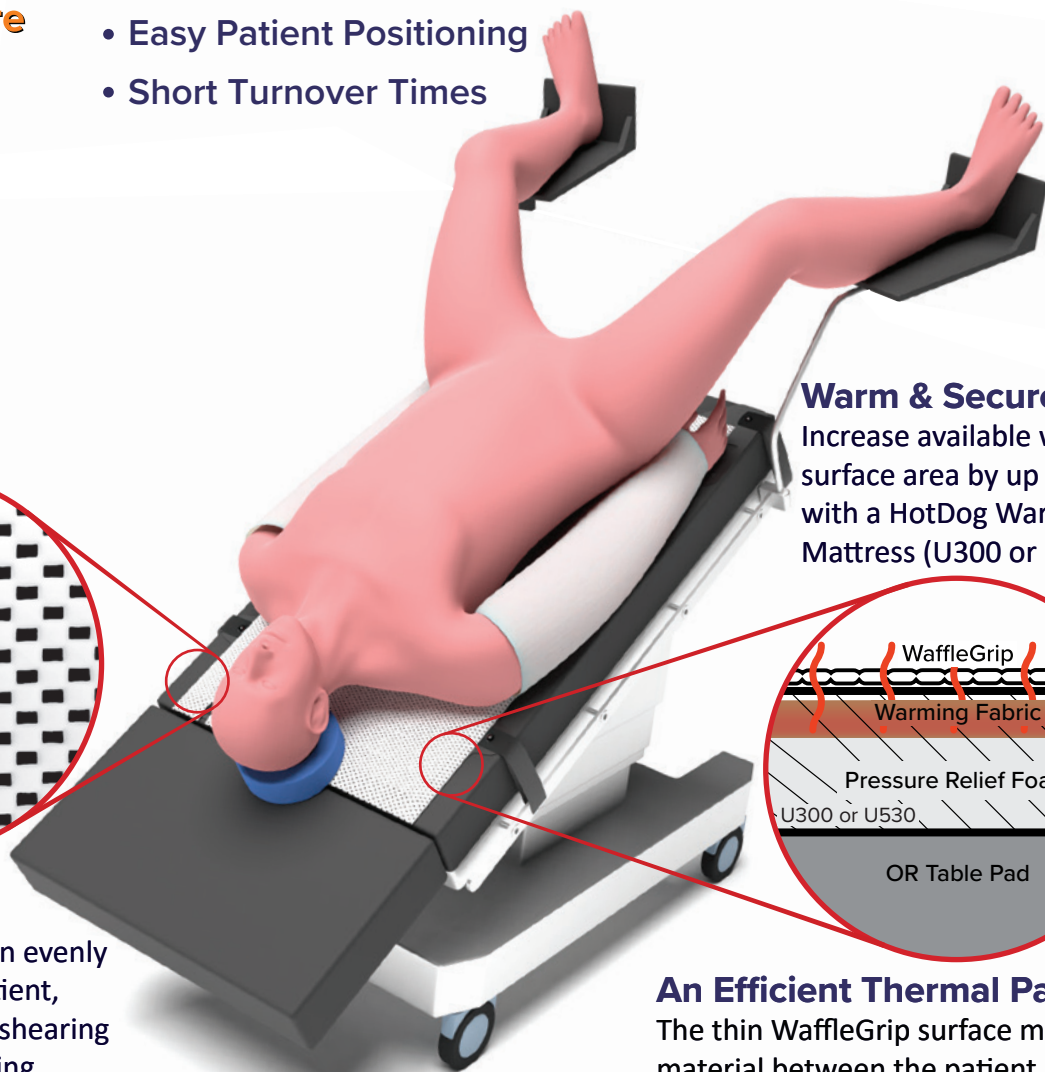

Warm & Secure

Increase available warming surface area by up to 3.5x with a HotDog Warming Mattress (U300 or U530).

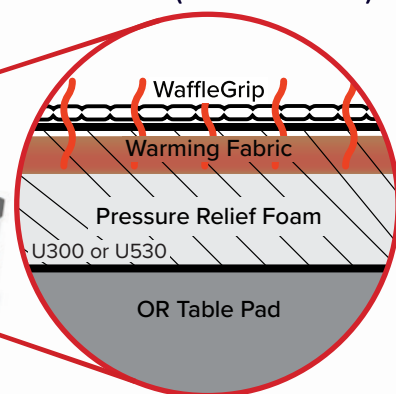

No Sliding

Unique waffle pattern evenly disperses grip on patient, decreasing localized shearing forces, while increasing mechanical hold and friction.

An Efficient Thermal Path

The thin WaffleGrip surface minimizes material between the patient and the HotDog warming fabric, supporting efficient conductive heat transfer.

Eleanor Markle, RN, Robotic Program Coordinator

“[WaffleGrip] is easy to use, easy to reposition, easy to educate staff and physicians on, and the warming is an additional benefit. Set-up is fast, so it doesn't delay OR turnover. It even comes with a slide sheet, making it very easy to reposition patients...Housekeeping loves, loves, loves it!”

Integration Optimizes Efficiency, Safety, & Outcomes



Efficient Underbody Warming

Trendelenburg positioning limits the surface area available for warming. WaffleGrip's thin profile enables underbody warming across up to 3.5x more surface area while supporting efficient heat transfer and secure patient contact.



Improve Securement

AORN guidelines recommend friction-based solutions for Steep T securement to reduce the risk of brachial plexus injury. WaffleGrip's anti-slip technology is suitable for high-BMI patients.



Streamlined Workflow

The integrated slide sheet allows effortless positioning to improve staff safety, and the perineal drape helps manage fluids...all while the patient remains Warm & Secure.




Cost Savings


The U300 and U530 mattresses are reusable. Only the disposable WaffleGrip kit is replaced after each procedure, while efficient setup and repositioning help save valuable staff time.

Easy Setup & Repositioning


Easily secure the U300 or U530 HotDog Mattress with intuitive, durable straps.

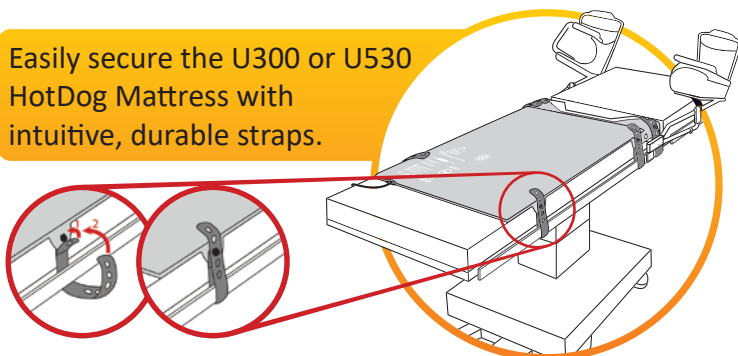


The BackSaver Slide Sheet allows you to easily slide the patient into position without straining your back. The Slide Sheet is easily removed prior to tilting the table.

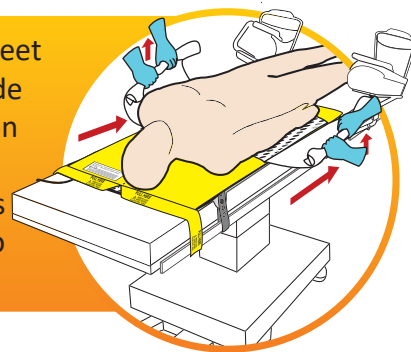

The Perineal Drape keeps the support surface clean by directing fluids to one designated area.

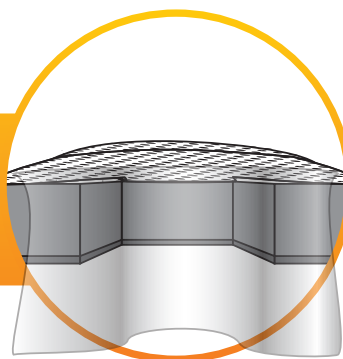

Keep the patient warm and securely positioned. (Shown with optional B500)

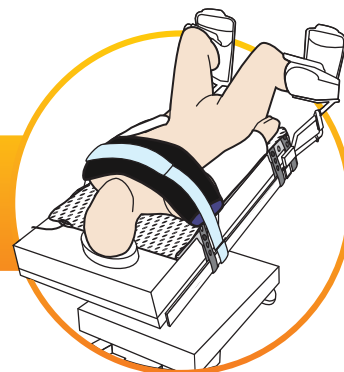

A301: For use with the U300 Warming Mattress

A302: For use with the U530 Warming Mattress + Return Electrode

Each package includes 10 complete disposable WaffleGrip Kits: a WaffleGrip pad with perineal drape, BackSaver Slide Sheet, and non-woven draw sheet with WaffleGrip.



**Augustine
Surgical Inc.**
Making Surgery Safer

Must be used in conjunction with HotDog Patient Warming System.

www.hotdogwarming.com/patents

©2026 Augustine Temperature Management. All Rights Reserved.

Rx Only

M214EN Rev E 07/26

6581 City West Parkway
Eden Prairie, MN 55344
952-465-3500 or 866-484-3505
email: info@aug Surg.com